

User's Manual

Aqua™ Shower Trolley BA M14 190 & BA M14 214



Contents:

1. Introduction
2. Area of Use
3. Taking Aqua into use
4. User Instructions
5. Maintenance
6. Trouble Shooting
7. Technical Data
8. Warranty

1. Introduction:

Aqua, electrically height-adjustable shower trolley

- 125mm castors with foot-operated central brake
- Folding guardrails in both sides
- Hand control for electrical functions
- Integrated battery and wall charging station
- Optional Equipment: Back-up battery

The Aqua Shower Trolley provides for problem free showering both for the carer and weak patients, who cannot sit. The trolley is comfortable thanks to the foldable guardrails, which are infinitely variable. This means that most patients obtain an optimal body support. The optional extra long trolley means that there is always ample access for head and foot care, even with very tall users.


The electrical lifting function provides optimal working positions for the carer both during showering and transfer to & from the trolley. The trolley can be tilted electrically. This increases comfort and safety. If needed a patient can quickly be put into either Trendelenburg or Anti-Trendelenburg position.

2. Area of Use: For washing / showering of weak patients.

3. Taking Aqua into use:

Upon reception: Unpack the trolley and check it thoroughly for transport damages. Put the plug into a wall socket and control that all electrical functions work. In case of damages or malfunctioning, contact your dealer **prior** to taking the trolley into use.

4. User Instructions:

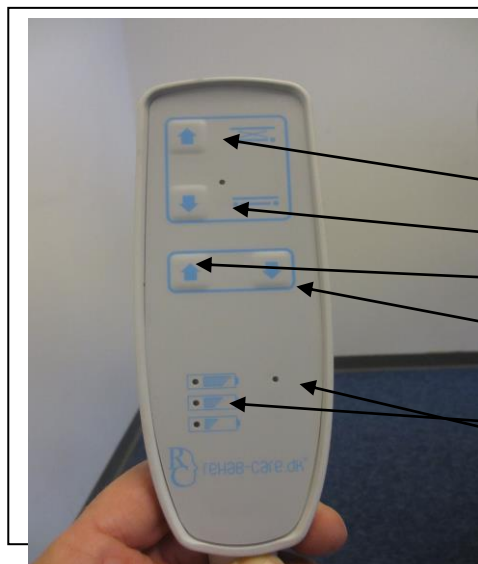


The electrical functions of the trolley are operated on the hand control (see illustration).

1. Raises the trolley
2. Lowers the trolley
3. Tilt the trolley in Trendelenburg
4. Tilt the trolley back to horizontal
5. Battery status
6. Indicate re-charge

All four driving wheels on the bearing are provided with central brake that can be operated with the foot.


Pressing the pedal with the red mark - brakes all 4 wheels.

Pressing the pedal with a green dot - activated directional lock on the 2 wheels in the foot.

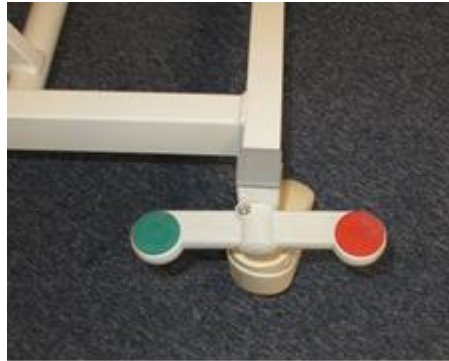

Adjusting the guardrails:

1. The guardrails can be adjusted when the release handle is pushed down.
2. Move rail to the desired position and let go.
3. Adjustment is step less.

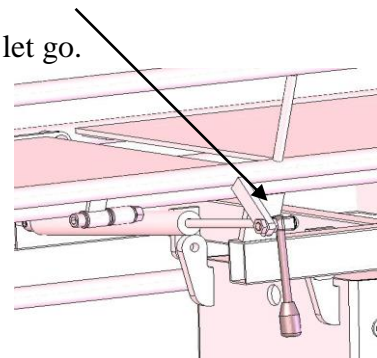

Head- & Footboards:

Can be removed for additional length.

The head- & footboards are fixated with thumbscrews underneath the mattress at both ends.



5. Maintenance:

Keep the trolley free of dirt, hair, soap etc.

Moving parts should regularly be controlled, tightened and lubricated with acid-free Vaseline.

Cleaning:

Clean the trolley with a soap water solution and a brush/sponge.

Wipe the trolley thoroughly dry after washing.

Cleaning agents:

-The trolley: Household detergents, **excluding alcohol & chloride agents.**

-Mattress: Only soap water

The trolley should never be cleaned in autoclaves or put into hospital washing streets!

A yearly Service inspection is recommended.

Contact your dealer for a service contract.

6. Trouble Shooting:

Malfunction of electrical functions:

1. Has the plug been inserted into the wall socket?
2. Does the wall socket supply power? Try another socket
3. Has the battery been charged?
4. Are all plugs driven all the way home in the control box? (situated under the undercarriage).
5. Is there any visual deformation to the trolley?
6. Are there any visual cracks in cables?

For further assistance, contact your dealer.

7. Technical Data:

Measurements & Weight:

• Length, external (Long version)	2430 mm
• Length, internal (Long version)	2140 mm
• Length, external (Short version)	2190 mm
• Length, internal (Short version)	1900 mm
• Width, internal	550 mm
• Total width, guardrails vertical	780 mm
• Total width, guardrails horizontal	1000 mm
• Total width, 1 guardrail horizontal	890 mm
• Height of guardrails	240 mm
• Length, head section	650 mm
• Castors	125 mm
• Height of mattress, trolley horizontal	430-960 mm
• Max height (Trendelenburg)	1340 - 1410 mm
• Minimum clearance	130 mm
• Length/width, undercarriage	1390/670 mm
• Thickness of mattress	30 mm
• Length of drain tube	1000 mm
• Weight	80kg
• Power needs	230V~50Hz.
• Max user weight	200kg

8. Warranty:

Warranty is granted on:

1. Frame & Mattress: 2 years
2. Other parts (castors, batteries, actuators etc.): 1 year

If the trolley is incorrectly used, the warranty is not valid.

Items numbers:

Aqua™ Shower Trolley, 190 cm length	BA M14 190
Aqua™ Shower Trolley, 214 cm length	BA M14 214